


DSC High Pressure Capsule Kit



Operator's Guide

Issued December 2010
Rev. K

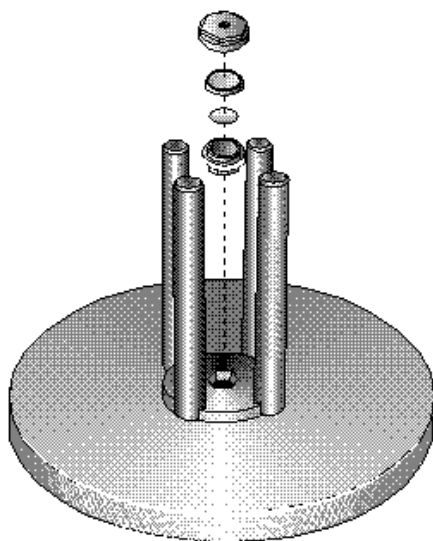

Notice

The material contained in this manual, and in the online help for the software used to support this instrument, is believed adequate for the intended use of the instrument. If the instrument or procedures are used for purposes other than those specified herein, confirmation of their suitability must be obtained from TA Instruments. Otherwise, TA Instruments does not guarantee any results and assumes no obligation or liability. TA Instruments also reserves the right to revise this document and to make changes without notice.

TA Instruments may have patents, patent applications, trademarks, copyrights, or other intellectual property covering subject matter in this document. Except as expressly provided in written license agreement from TA Instrument, the furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property.

TA Instruments Operating Software, as well as Module, Data Analysis, and Utility Software and their associated manuals and online help, are proprietary and copyrighted by TA Instruments. Purchasers are granted a license to use these software programs on the module and controller with which they were purchased. These programs may not be duplicated by the purchaser without the prior written consent of TA Instruments. Each licensed program shall remain the exclusive property of TA Instruments, and no rights or licenses are granted to the purchaser other than as specified above.

Important: TA Instruments Manual Supplement

Click on the links below to access important information supplemental to this Getting Started Guide.

- [TA Instruments Trademarks](#)
- [TA Instruments Patents](#)
- [Other Trademarks](#)
- [TA Instruments End-User License Agreement](#)
- [TA Instruments Offices](#)

Notice

Please read carefully:

Except for any express warranties contained herein, TA Instruments disclaims all other warranties, including, without limitation, all implied warranties of merchantability and fitness for a particular purpose.

TA Instruments warrants that, as delivered to you, the product conforms to our quality standards for use in the TA Instruments Differential Scanning Calorimeters. Your exclusive remedy for TA Instruments' breach of any warranty will be, at TA Instruments' sole option, to use reasonable commercial efforts to attempt to correct the problem, to replace the defective product, or to refund the purchase price and terminate this agreement.

This agreement sets forth the entire agreement and understanding of the parties with regard to this transaction and supersedes all prior agreements, discussions, and understanding between them, whether oral or written.



WARNING: These capsules should not be used for analysis of thermally unstable or explosive materials. You must always use extreme care for your own and others' safety when handling materials that can decompose violently.

Table of Contents

Important: TA Instruments Manual Supplement	3
Notice	4
Please read carefully:	4
Table of Contents	5
Overview	7
Safety	7
Safety Label	7
Cleaning the Capsules	9
Summary of High Pressure Capsule Safety Devices	10
Preparing a Capsule	11
Sealing and Weighing the Sample	11
Opening the Capsule	13
Calibrating the DSC	14
DSC Q200/Q100/Q20A or Q2000/Q1000 Instruments	14
DSC Q20/Q10 Instrument	14
DSC 2920/2910/2010 Instruments	16
Running Experiments	17
Maintaining the Capsules and Sealing Tool	18
Stainless Steel Capsules	18
Sealing Tool Maintenance	18
Specifications	18
Replacement Parts	19
Index	21

Overview

This booklet describes the procedures needed to prepare high pressure sample capsules for use in the DSC cell. Be sure to read the following section on safety carefully before proceeding further.

Safety

The high pressure capsules have been evaluated as to European Union requirements in EN61010-1/1993 + A2/1995 and EN61010-2-010/1994 and have been found to fall below the pressure x volume product that is used to determine applicability.

Safety Label



The label shown at the left is displayed on the metal bell jar that comes with the DSC High-Pressure Capsule Kit for your protection. This label indicates that a hot surface may be present. Take care not to touch this area or to allow any material that may melt or burn to come in contact with this hot surface.

The following notes, cautions, and warnings are provided to prevent accident and ensure safe laboratory practice. Read this section, beginning on the next page, carefully before using the materials provided in this kit, and follow all instructions. For protection, wear safety glasses at all times.



WARNING: High pressure capsules are designed for evaluation of materials *under self-generated atmospheres at pressures up to 10 MPa (1450 psi)*. As such, they are *not* intended to be a replacement for the pressure DSC cell, which is designed to evaluate materials under high pressures of specific purge gas (usually inert or reactive). Rather, the capsules are primarily intended for the evaluation of materials for which suppression of volatilization of water or solvent (or sublimation) is required to obtain good heat flow results.

With a 10- μ L (nominal) aqueous sample,* these capsules can be used to about 300 °C before the upper limit is exceeded. However, since other samples may build up pressure more rapidly on heating and, more important, since the final rupture (failure) of the capsules at the upper pressure limit could result in damage to the DSC cell, it is strongly recommended that you use caution in deciding whether or not to evaluate a material in the high pressure capsules. Highly energetic materials such as pyrotechnics should never be run in these capsules.

*Do not exceed sample volume capacity (35 μ L) of the capsule (*i.e.*, do not fill capsule completely with sample).



WARNING: Do not touch the inside of the cell, or lean over it when inserting or removing a capsule. Handle cooled capsules with care. Gases produced during decomposition reactions at high temperatures may not condense when cooled; therefore, the capsule may remain under pressure at ambient temperature.



WARNING: Properly sealed capsules may release pressure rapidly when internal pressure exceeds capsule-seal pressure capability. These capsules should not be used for analysis of thermally unstable or explosive materials.

DSC Q Series or Discovery DSC™: The high pressure capsule can be used with the Discovery DSC, Q2000/Q1000, or Q200/Q100/Q20A models configured with a Finned Air Cooling System (FACS), RCS, or LNCS. The AutoLid must be in position during experiments. For the DSC Q20/Q10 instrument, the manual lid and safety clamp must be in position during experiments.

DSC 2920/2910/2010: Be sure to use the supplied safety devices (metal bell jar for cells without a heat exchanger, or safety lid for cells using a heat exchanger) when running any experiments with the high pressure capsules.



WARNING: The safety lid provided will not work for 2920/2910/2010 Series RCS (Refrigerated Cooling System) heat exchangers that have been modified with a “prototype” heat transfer sleeve. RCS units with serial numbers lower than 1641 may contain this prototype sleeve. Contact our service department at (302) 427-4050 or your local TA Instruments Representative for further details on how to modify those RCS units for use with high pressure capsules. There are no concerns regarding the use of the LNCA (Liquid Nitrogen Cooling Accessory) heat exchanger with the DSC high pressure capsules.



WARNING: The DSC high pressure capsules should not be used with an older-style Mechanical Cooling Accessory (MCA). The MCA (PN 990460.901/.902), which was discontinued as a product in 1995, does not work with the safety lid included in the DSC High Pressure Capsule Kit.



CAUTION: Properly sealed capsules can withstand internal pressures up to 10 MPa (1450 psi), with temperatures up to 300 °C. Improperly sealed capsules can leak before reaching this pressure.



CAUTION: Heating rates should be 10 °C/minute or less because of the high heat capacity of the capsule. Keep the sample mass as low as practical.

Cleaning the Capsules



CAUTION: Handle the capsule lid, capsule bottom, seal, sample, and assembled caps, covers, and sealing disks with tweezers or other suitable tool. The body oils from your hands can give erroneous data.

Like the other pans provided for use with TA Instruments DSC systems, the high pressure capsules are manufactured to high quality standards, including cleaning to remove contaminants that might be present from the manufacturing process. For most applications, these pans can be used as received; however, if the pans are used for high sensitivity experiments, an additional cleaning process is recommended before use. This procedure is taken from Appendix A of ASTM standard E1858, *Test Method for Oxidative Induction Time of Hydrocarbons by Differential Scanning Calorimeters*.

Follow the procedure given here to clean TA Instruments DSC high pressure capsules:

1. Place a maximum of 20 capsules in a 250-mL Erlenmeyer flask that has been fitted with a glass stopper.
2. Add approximately 150 mL of reagent grade toluene (enough to cover the capsules).
3. Swirl the flask, containing the capsules and toluene, for 0.5 to 2.0 minutes.
4. Let the flask stand for at least 1 minute.
5. Decant the toluene out of the flask.
6. Repeat steps 1 through 5.
7. Add approximately 150 mL of reagent grade acetone after the second toluene wash.
8. Swirl the flask, containing the pans and acetone, for 0.5 to 2.0 minutes.
9. Let the flask stand for at least 1 minute.
10. Decant the acetone out of the flask.
11. Repeat steps 7 through 10.
12. Rotate the flask—so that no capsules adhere to the bottom or sides of the flask—as you flow nitrogen at 150 to 200 mL/min over the wet capsules to drive off the excess solvent. This should take approximately 5 to 6 minutes.
13. Return the cleaned capsules to their storage container, and record the date they were cleaned.

Summary of High Pressure Capsule Safety Devices

<i>Instrument</i>	<i>High Pressure Capsule Safety Device</i>
Discovery DSC, Q2000/Q1000, Q200/Q100/Q20A	AutoLid (See page 14.)
DSC Q20/Q10	Manual Lid with Safety Clamp (See pages 14 and 15.)
DSC 2920/2910/2010 —without heat exchanger	Metal Bell Jar (See page 16.)
DSC 2920/2910/2010 —with heat exchanger —with RCS units having serial numbers above 1641	Safety Lid (See page 16.)
DSC 2920/2910/2010 —with MCA units having PN 990460.901/.902 (discontinued in 1995) —with RCS units having serial numbers below 1641	CANNOT USE HIGH PRESSURE CAPSULES

Preparing a Capsule



CAUTION: Handle the capsule lid, capsule bottom, seal, sample, and assembled caps, covers, and sealing disks with tweezers or other suitable tool. The body oils from your hands can give erroneous data.

Refer to Figure 1 and follow the instructions in this section to prepare sample capsules for the DSC high pressure capsule experiments.

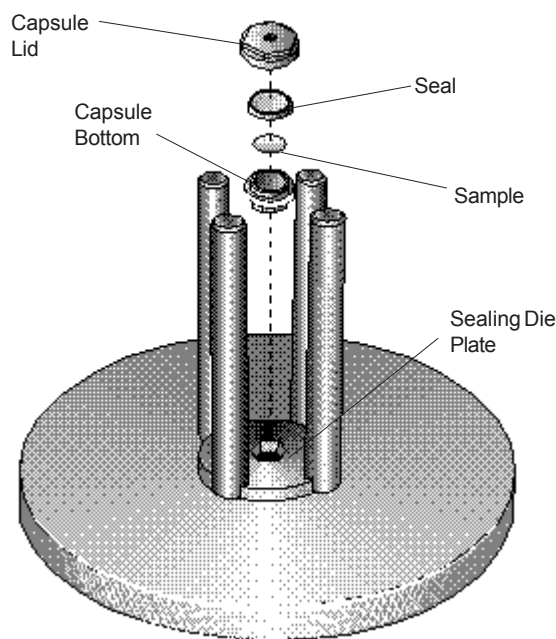


Figure 1
Parts of DSC High Pressure Capsule Sealing Assembly

Sealing and Weighing the Sample

Before you can enter the sample size in your experimental parameters information, you will need to weigh the sample.

NOTE: Try to keep the sample mass as low as practical..

To weigh samples, follow these directions:

1. Obtain the total weight of the capsule bottom, capsule lid, and seal.
2. Using tweezers, place a capsule bottom in the hexagonal hole of the sealing plate (see Figure 1).
3. Place the sample carefully in the capsule bottom; then place the seal, cup downward, on the capsule (see Figure 1).

4. Using tweezers, place the capsule lid over the seal.
5. Align the sealing tool over the hexagonal lid as shown in Figure 2 here.
6. Turn the sealing tool clockwise, while applying light pressure, until the sealing tool slips—you will hear a click.
7. Remove the sealing tool from the sealing die plate, and take out the sealed capsule.

NOTE: When handling the sealed capsule, be sure to hold it level, especially when using liquid and viscous samples. The sample will then remain in the pan, and will not adhere to the underside of the lid.

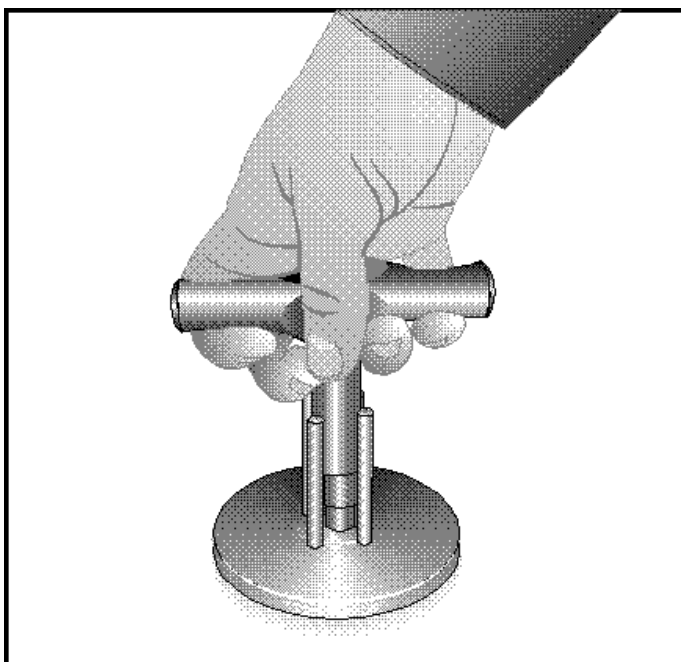


Figure 2
Using the Sealing Tool

8. Weigh the capsule with the seal and sample. Determine the sample weight by subtracting the weight of the empty capsule (found in step 1) from the total weight:

$$\text{Weight}_{\text{step 8}} - \text{Weight}_{\text{step 1}} = \text{Sample weight}$$

The high pressure capsule is now ready to run in your DSC cell.

Opening the Capsule

To open the capsule, reverse the procedures listed previously for sealing. This allows you to remove the sample and reuse the capsule. Seals must be discarded after use.



WARNING: High pressure capsules may retain pressure after material decomposition. Make sure the capsule is at room temperature before opening it. Open the capsule slowly and carefully to allow the seal to open and vent the pressure.

Calibrating the DSC

In order to get accurate experimental results, your DSC system should be recalibrated for cell constant and temperature using the sealed high pressure capsules. Follow these directions for calibration. Refer to the appropriate section for your instrument.

Discovery DSC Q200/Q100/Q20A or Q2000/Q1000 Instruments

1. Use an empty sealed capsule as a reference. Place it in the cell.
2. Weigh an 8 to 10 mg sample of indium.
3. Place a thin layer of alumina (about 20 mg) in the sample capsule. Then place the indium sample on top of the alumina.

NOTE: Normally DSC heat flow and temperature calibration are accomplished by evaluating the calibration material (usually indium) under the same experimental conditions as the subsequent sample materials. With the increased resolution of the Q DSC T4 heat flow, however, the calibration conditions must be changed slightly to account for the small indium sample mass relative to the large mass of the low thermal conductivity stainless steel high pressure capsule. The addition of alumina helps counter balance that difference and provides a calibration better suited to the kinds of material run in these capsules.

No alumina is necessary for T1 heat flow calibration.

4. Seal the capsule and place it inside the cell.
5. Close the AutoLid before proceeding. This step is very important for safety reasons.
6. Use the recommended calibration methods to calibrate the system (refer to the DSC and the instrument control online documentation for information). Make sure that calibration of the Tzero cell resistance and capacitance values are performed.

NOTE: The Calibration Wizard can be used to conveniently calibrate the DSC Q Series instruments.

7. Return the instrument to standard mode before running your experiments.

DSC Q20/Q10 Instrument

1. Use an empty sealed capsule as a reference. Place it in the cell.
2. Weigh an 8 to 10 mg sample of indium.
3. Seal the capsule and place it inside the cell.
4. Install the safety clamp on the manual lid as follows:
 - a. Unscrew the knob from the top of the manual lid.
 - b. Place the hole in the middle of the safety clamp over the exposed screw.



- c. Screw the knob back onto the manual lid and tighten.
5. Squeeze the tabs on the safety clamp toward the center (as shown in the figure here) while you lower the manual lid onto the cell.

IMPORTANT: Make sure the hooks at either end are fully engaged in the vent holes as shown. Pull up slightly on the lid to make sure it is secure.



WARNING: Do not attempt to run the DSC Q20/Q10 using the high pressure capsules unless the safety clamp is securely in place as shown.

6. Use the recommended calibration methods to calibrate the system. The Calibration Wizard can be used to conveniently calibrate the DSC Q20/Q10. Only T1 heat flow calibration is possible on the Q20/Q10.
7. Return the instrument to standard mode.

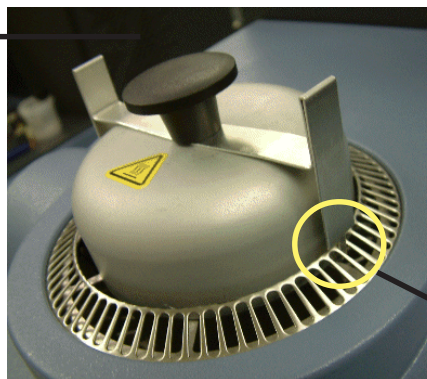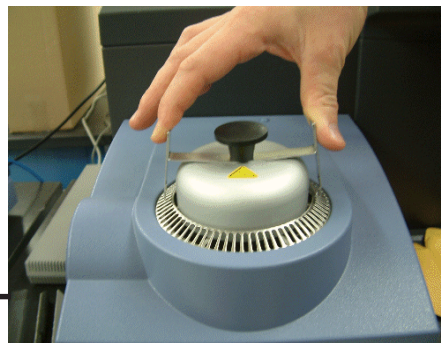

Position
Hooks
Securely

DSC 2920/2910/2010 Instruments

1. Use an empty sealed capsule as a reference. Place it in the cell.
2. Seal an 8 to 10 mg indium metal sample in a high pressure capsule and place it in the cell.
3. Install the *metal bell jar* for cell without a heat exchanger. See Figure 3. Install the *safety lid* for cell using a heat

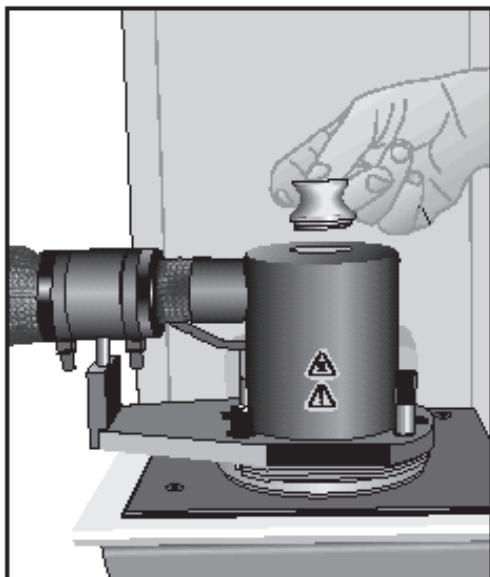


Figure 4
Using the Safety Lid with the 2920/2910/
2010 DSC Cell
and the Heat Exchanger

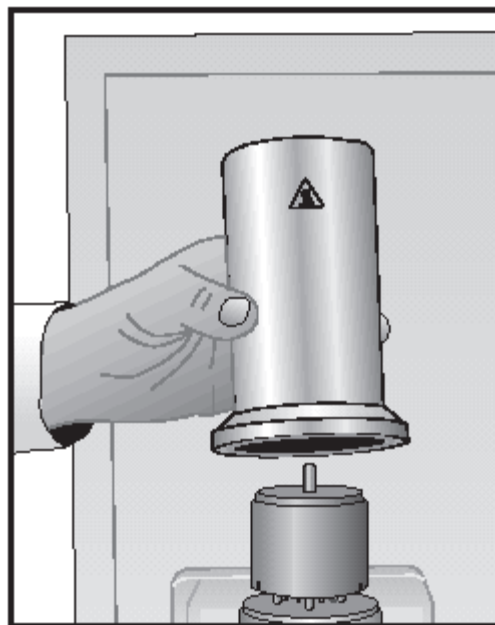


Figure 3
Using the Metal Bell Jar with the 2920/
2910/2010 DSC Cell

exchanger. See Figure 4.



WARNING: Properly sealed capsules may release pressure rapidly when internal pressure exceeds capsule-seal pressure capability. Be sure to use the supplied safety devices (metal bell jar for cells without a heat exchanger, or safety lid for cells using a heat exchanger) when running any experiments or calibrating with the high pressure capsules.

Make sure that you read the Warnings beginning on page 7 BEFORE proceeding.

4. Use standard calibration methods to calibrate the system. If additional temperature calibration points are required, make sure that you use the high pressure capsule when gathering data.
5. Return the instrument to standard mode before running your experiments.

Running Experiments

After you calibrate the DSC using the high pressure capsule and return the instrument to standard mode, load the sample and reference capsules into the DSC cell.



WARNING: Properly sealed capsules may release pressure rapidly when internal pressure exceeds capsule-seal pressure capability.

DSC 2920/2910/2010 Instruments: Be sure to use the supplied safety devices, as shown in Figures 3 and 4 on the previous page, when running any experiments or calibrating with the high pressure capsules.

Discovery DSC, Q200/Q100/Q20A, and Q2000/Q1000 Instruments: Be sure to close the AutoLid before beginning the experiments using high pressure capsules.

DSC Q20/Q10 Instruments: Make sure the safety clamp is securely installed on the manual lid before beginning the experiments using high pressure capsules.



WARNING: Do not use the high pressure capsules in a DSC 2920 dual sample cell. The high pressure capsules are not designed to be used with that cell.



WARNING: Make sure that you read the Warnings beginning on page 7 BEFORE proceeding.

Maintaining the Capsules and Sealing Tool

The capsule bottom and lid are normally reusable after being run on the instrument. Open the capsule as directed, remove the sample, and clean the capsule after each use. The seal must be replaced after each use of the capsule.

Stainless Steel Capsules

The stainless steel capsules can react differently from other types of capsules; they:

- May discolor and distort at high temperatures.
- May be cleaned with a solvent that does not affect the steel.

Sealing Tool Maintenance

Wipe the sealing tool and sealing die plate clean with a soft cloth that has been dampened with a dilute laboratory detergent solution.

The hexagonal hole in the sealing die plate should be periodically checked for debris from the capsules. Use a soft tissue or cloth to clean the sealing die plate as needed.

Specifications

Table 1 contains the technical specifications for High Pressure Capsules and seals.

Table 1
High Pressure Capsule Specifications

Pressure capability	10 MPa (1450 psi)
Temperature limit for aqueous solutions	300 °C
Capsule Volume	35 μ L (max)
Material	440A SST
Seal material	Gold-plated copper

Replacement Parts

Table 2 contains the replacement parts for High Pressure Capsules and seals.

Table2
List of DSC High Pressure Capsule Parts

<i>Part Number</i>	<i>Description</i>
900814.901	DSC High Pressure Capsule Seals (qty 20)
900815.901	DSC High Pressure Capsules (qty 5)

Index

A

alumina 17

C

calibrating DSC
for high pressure capsules 17
DSC Q200/Q100 or Q2000/Q1000 instru-
ments 17

Calibration Wizard 17

capsule
volume capacity 8

capsules
calibrating DSC 17
cleaning 22
after use 22
before use 11
handling precautions 11
handling when sealed 15
heat capacity 9
improperly sealed 9
internal pressure 9
maintenance 22
material 23
opening 16
parts 24
preparing 13
releasing pressure 20
reordering 24
reusing 22
sealing 14
specifications 23, 24
stainless steel 22
temperature limits 8
volume 23
weighing 15

cleaning
sealing die plate 22

D

decompose
materials 5

decomposition reactions 8

DSC instrument
calibrating for high pressure capsules 17

DSC Q20/Q10 safety clamp 18

DSC Q200/Q100 and Q2000/Q1000
high pressure capsules 17

E

European Union
safety requirements 7

experiments with high pressure capsules 20

G

gases 8

H

heat exchanger
safety lid for 21

heating rate limits 9

hexagonal lid 14

high pressure capsule kit
safety devices 18, 20

M

materials
decomposing 5
self-generated atmosphere 7

metal bell jar 21
safety label 7

N

notice of warranty 5

P

parts
high pressure capsule 24

pressure
internal 8

pressure capability 23

pressure limits 7

R

RCS (Refrigerated Cooling System)
safety lid 9

reference 17, 18, 19

running experiments with high pressure capsules 5

S

safety 7
CE compliance 7
heating rate limits 9
label 7
pressure/temperature limits 7

safety clamp 18
installation 18

safety lid
for DSC cells with heat exchanger 21

sample
decomposition 16
volume capacity 8
weighing 14

seal(s)
material 23

sealing die plate
maintenance 22

sealing tool
maintenance 22

sealing tool/assembly
cleaning 22
maintenance 22
parts 13
using 14

seals
reordering 24

T

T1 heat flow 18

T4 heat flow 17

temperature limit 23

V

volume
capacity of capsule 8

W

warranties 5